

# waterpik® waterflosser®



## QUICK START GUIDE

Review full instruction book prior to first use

1

### FIRST USE

Fill reservoir with warm water and replace it firmly on the base.

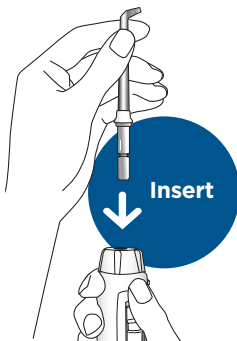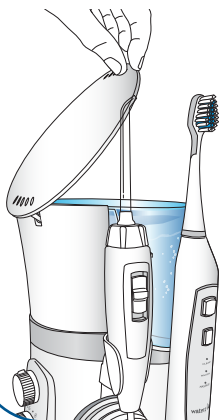

Insert tip firmly until it clicks into place.

Turn pressure to high, point tip into sink, turn power and water switches on until water flows to prime unit. Then turn off.



## QUICK START GUIDE

Review full instruction book prior to first use

### 2 START WATER FLOSSING

**TO AVOID MESS: PLACE THE TIP IN YOUR MOUTH BEFORE TURNING UNIT ON**


Turn power and water on.  
Start on low pressure and increase to your liking.

Water ON



### 3 LEAN OVER SINK

Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

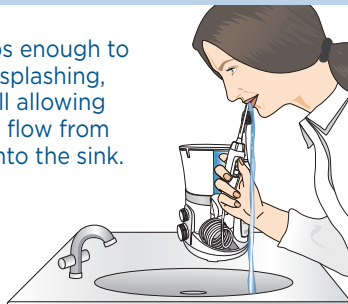

### 4 TECHNIQUE

Aim water at the gumline at a 90° angle. Follow the gumline and pause briefly between teeth.

